



Electrophorus


Handle & Plates

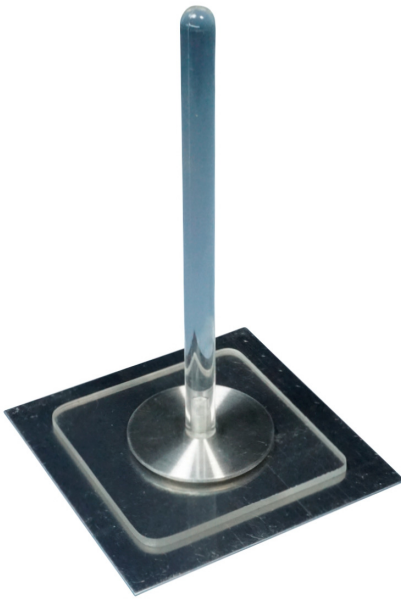

EM1780-001
(50mm Diameter)

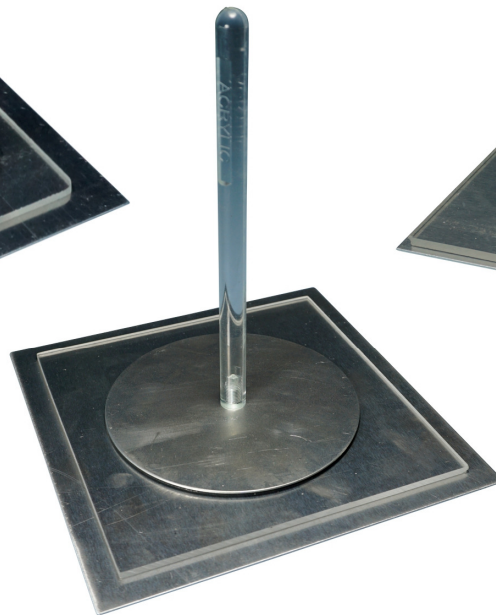

EM1781-001
(75mm Diameter)

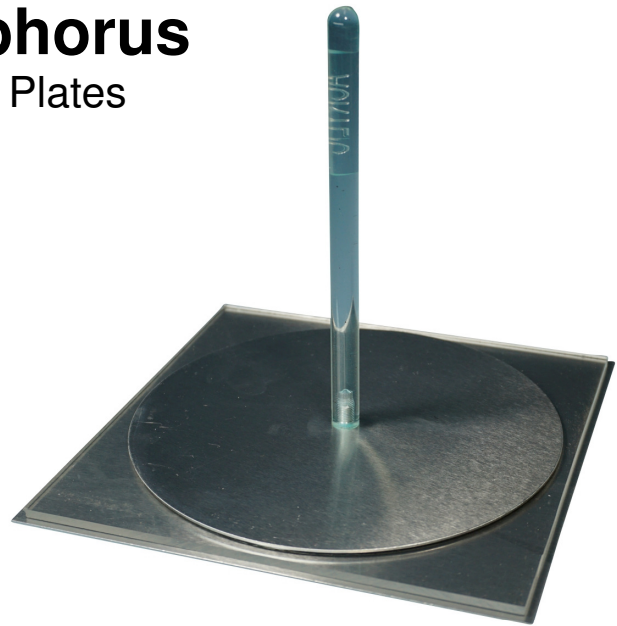

EM1782-001
(150mm Diameter)

Description:

The classic Electrophorus is a metal disc with insulating handle, a metal base plate and an insulating plate that rests on the base plate. When an electric charge is applied between the base plate and the disc, a capacitor is formed and the plate can be removed with an electric charge residing.

| Code: | Disc Diameter: | Weight: |
|------------|----------------|---------|
| EM1780-001 | 50mm | 100g |
| EM1781-001 | 75mm | 150g |
| EM1782-001 | 150mm | 200g |

Designed and manufactured in Australia

INDUSTRIAL EQUIPMENT & CONTROL PTY.LTD.

61-65 McClure St. Thornbury. 3071 Melbourne. Australia

Tel: 61 (0)3 9497 2555

Email: iec@iecpl.com.au

www.iecpl.com.au